

Svea TW



Indoor Use
Shock Absorber
Easy to Adjust



The Svea TW is recommended for indoor use. This sturdy and safe trolley walker features a hydraulic gas spring that acts as a shock absorber while walking, slowing down the vibrations caused by movement and ensuring that the wheels remain in permanent contact with the ground.

Hydraulic Shock Absorber

The Svea TW features a unique frame design that incorporates a hydraulic gas spring to lessen the vibrations caused by movement while walking. This is especially helpful to users that suffer from pain or neurological disorders. The wide, firm table offers the extra support required and an optional stabilizer wheel can be added for more support.

Unique Handles

The Svea TW features unique front-facing, angle adjustable handles that allow users to find the optimal handle position for comfort and stability. Adjustment is easy with integrated measurements right in the handle tube. The ergonomic design of the grips

ensures even distribution of weight. The Staffan TW can also be equipped with a one-hand brake to allow a user to control hand functions on either the left or right hand side.

High Capacity

The Svea TW is available in two widths and two heights and combines maximum weight bearing relief with the manoeuvrability of a regular rollator. With a high weight capacity of 150kg/330.7lb, the Svea TW is easy to adjust and use in most home environments.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories/Features



Stabilizer Wheel

For added stability while walking. Part no. 50609-Svea



Security Lock

Flex cord that protects rollator from theft while unattended. Part no. 91092




One-Hand Brake

Allows user to control hand functions with one hand (installed left or right). Part no. 10653



Cane holder

Simplifies transport of a cane or umbrella. Cannot be used when one-hand brake is installed. Part no. 10603

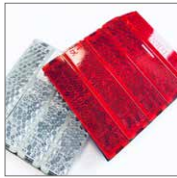

Reflector Kit

Includes 30 self-adhesive reflectors. Part no. 91090



Hydraulic Gas Spring

Acts as a shock absorber while walking. Part no. 50552 (cylinder)



Light kit

Illuminate your way while walking. Part no. 91091

Note: Not all accessories are available for sale in all countries.

Technical Specifications

	Svea TW Low 59	Svea TW Nor 59	Svea TW Low 68	Svea TW Nor 68
Part no.	53211	53221	53411	53421
Handle Height (Min - Max)	82cm - 107cm / 32.3" - 42.1"	101cm - 126cm / 39.7" - 49.6"	82cm - 107cm / 32.3" - 42.1"	101cm - 126cm / 39.7" - 49.6"
Overall Width	59cm / 23.2"	59cm / 23.2"	68cm / 26.8"	68cm / 26.8"
Wheel Size	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"
Product Weight	14.9kg / 32.8lb	15.1kg / 33.3lb	15.1kg / 33.3lb	15.3kg / 33.7lb
Maximum User Weight	150kg/330.7lb	150kg/330.7lb	150kg/330.7lb	150kg/330.7lb

Follow Us on Social Media